SURGICAL INSTRUMENTS

TYPES OF SURGICAL PACKS

× OPHTHALMIC PACK (Eyes)

× LACERATION PACKS (Small)

× SURGICAL PACKS (Large)

ALL PACKS SHOULD BE

- Double Wrapped(1 Cloth Drape Inside, 1 Paper Drape Outside)
- Include 30 pieces of Gauze
- × Include a Sterile Strip Indicator



INSTRUMENTS INCLUDED IN OPHTHALMIC PACK

- × 4 Towel Clamps
- × 1Thumb Forcep
- × 1 Blade Handle
- × 1 Needle Driver
- × 1 Castroviejo
- × 1 Mayo Scissors
- × 1MetzenbaumScissors

- × 2 Pairs Iris Scissors
- × 1Blepharstats
- × 2 Curved Mosquitos
- × 2 Straight Mosquitos
- × 3 Fixation Forceps
- × 24X24 Drape
- × 30 pc. Gauze
- × Indicator Strip



INSTRUMENTS INCLUDED IN LACERATION PACK

- × 4 Towel Clamps
- × 1Thumb Forcep
- × 1 Blade Handle
- × 1 Needle Holder
- × 1 Mayo Scissor
- × 1Metzenbaum Scissor
- × 1 Sharp/Blunt Scissor

- × 3 Curved Mosquitos
- × 30 x 30 Drape
- × 30 pc. Gauze
- × OK Strip Indicator



INSTRUMENTS INCLUDED IN SURGICAL PACK

- × 6 Towel Clamps
- × 1Thumb Forcep
- × 1 Dressing Forcep
- × 1 Blade Handle
- × 1 Needle Driver
- × 1 Suction Tip
- × 1 Mayo Scissor
- × 30 pc Gauze & Ok Strip

- × 1Metzenbaum Scissor
- × 1 Sharp/Blunt Scissor
- × 3 Lg Curved Forceps
- × 1 Straight Med Forcep
- × 3 Med Curved Forceps
- × 1Straight Sm Forcep
- × 3 Sm Mosquitos
- × 2 Allis Forceps

WHAT IS A BUNT CLIP

Holds instruments together to keep them organized for the surgeon








NEEDLE DRIVER OR NEEDLE HOLDER



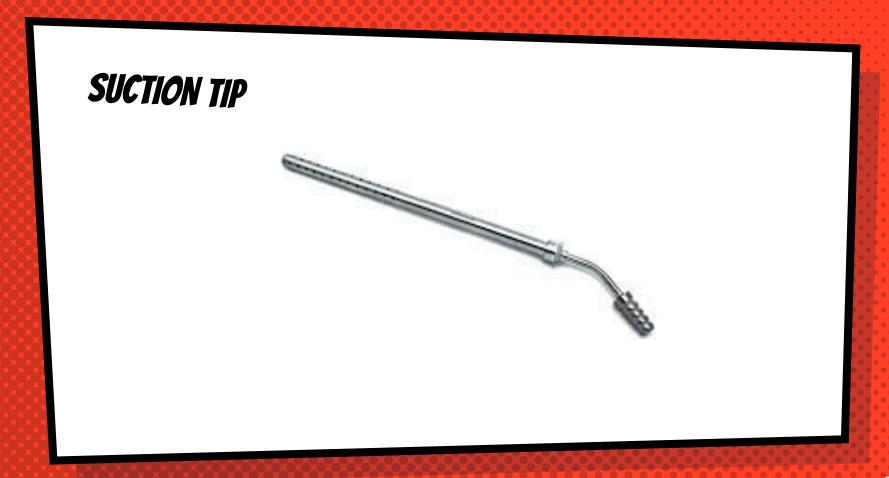

MAYO SCISSORS



SHARP/BLUNT SCISSORS



BACKHAUS TOWEL CLAMPS



ALLIS (THE ALLIGATOR) FORCEPS



KELLY FORCEPS (CURVED & STRAIGHT)

MOSQUITO FORCEPS (CURVED & STRAIGHT)



3X3 GAUZE, OK STRIP & PAPER DRAPE








LAP SPONGES & SUCTION TUBING





BLEPHARSTATS

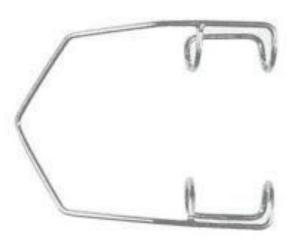

CASTROVIEJO NEEDLE HOLDERS



IRIS SCISSORS



BALFOUR RETRACTOR



GELPI RETRACTOR








GAS STERILIZATION

- Uses ethylene oxide to sterilize
- Must be certified to use it
- It is a colorless, odorless and POISONOUS gas that attacks the cellular proteins and nucleic acids of microorganisms.
- Used for plastic and rubber instruments

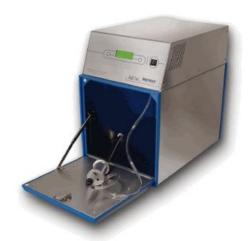

STEAM STERILIZATION

- Uses high temp & high pressure steam to sterilize
- Must use distilled H2O
- This is for packs & pouches
- Used for non-plastic instruments



OUR PROCESS IS EASY STERILIZE WRAP **CLEAN**

THAS!

Any questions?
You can find me at sgarcia@aavec.com

